

Medtronic

53401

Single Chamber Temporary External Pacemaker

Electromagnetic Compatibility Declaration

The following list includes trademarks or registered trademarks of Medtronic in the United States and possibly in other countries. All other trademarks are the property of their respective owners.

Medtronic

1 Electromagnetic compatibility declaration

The Medtronic Model 53401 Single Chamber Temporary External Pacemaker (hereafter referred to as the “temporary pacemaker”) is tested with the Medtronic 5433A and 5433V Patient Cables, and meets the requirements of IEC 60601-1-2:2007/AC:2010.

- Use of accessories and cables other than those that are specifically listed may result in increased emissions or decreased immunity of the temporary pacemaker.
- The temporary pacemaker should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, the temporary pacemaker should be observed to verify normal operation in the configuration in which it will be used.

Guidance and manufacturer’s declaration—electromagnetic emissions		
The temporary pacemaker is intended for use in the electromagnetic environment specified below. The customer or the user of the temporary pacemaker should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment—guidance
RF emissions CISPR 11	Group 1	The temporary pacemaker uses RF energy primarily for its internal function. Therefore, its RF emissions are very low and are not likely to cause interference in nearby electronic equipment. The temporary pacemaker is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class A	
Harmonic emissions IEC 61000-3-2	Not Applicable	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Not Applicable	

Guidance and manufacturer’s declaration—electromagnetic immunity			
The temporary pacemaker is intended for use in the electromagnetic environment specified below. The customer or the user of the temporary pacemaker should assure that it is used in such an environment.			
Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment—guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	Not applicable	Not applicable	
Surge IEC 61000-4-5	Not applicable	Not applicable	
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Not applicable	Not applicable	

Guidance and manufacturer's declaration—electromagnetic immunity			
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

Guidance and manufacturer's declaration—electromagnetic immunity			
The temporary pacemaker is intended for use in the electromagnetic environment specified below. The customer or the user of the temporary pacemaker should assure that it is used in such an environment.			
Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment—guidance
Conducted RF	3 V _{RMS} 150 kHz to 80 MHz outside ISM bands ^a	10 V	Portable and mobile RF communications equipment should be used no closer to any part of the temporary pacemaker, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 0.35 \sqrt{P}$
IEC 61000-4-6	10 V _{RMS} V _{RMS} MS 150 kHz to 80 MHz outside ISM bands ^a	10 V	$d = 1.2 \sqrt{P}$
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.5 GHz	10 V/m	$d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). ^b Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^c should be less than the compliance level in each frequency range. ^d Interference may occur in the vicinity of equipment marked with the following symbol: 

Guidance and manufacturer's declaration—electromagnetic immunity

Note 1: At 80 MHz and 800 MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

- ^a The ISM (industrial, scientific, and medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.
- ^b The compliance levels in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2.5 GHz are intended to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas. For this reason, an additional factor of 10/3 has been incorporated into the formulae used in calculating the recommended separation distance for transmitters in these frequency ranges.
- ^c Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the temporary pacemaker is used exceeds the applicable RF compliance level above, the temporary pacemaker should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the temporary pacemaker.
- ^d Over the frequency range of 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Recommended separation distances between portable and mobile RF communications equipment and the temporary pacemaker

The temporary pacemaker is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the temporary pacemaker can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the temporary pacemaker as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter	Separation distance according to frequency of transmitter			
	150 kHz to 80 MHz outside ISM bands $d = 0.35 \sqrt{P}$	150 kHz to 80 MHz in ISM bands $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$
0.01 W	0.04 m	0.12 m	0.12 m	0.23 m
0.1 W	0.11 m	0.38 m	0.38 m	0.73 m
1 W	0.35 m	1.2 m	1.2 m	2.3 m
10 W	1.11 m	3.79 m	3.79 m	7.27 m
100 W	3.5 m	12 m	12 m	23 m

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2: The ISM (industrial, scientific, and medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.

Recommended separation distances between portable and mobile RF communications equipment and the temporary pacemaker


Note 3: An additional factor of 10/3 has been incorporated into the formulae used in calculating the recommended separation distance for transmitters in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2.5 GHz to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas.

Note 4: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Medtronic



Medtronic, Inc.
710 Medtronic Parkway
Minneapolis, MN 55432
USA
www.medtronic.com
+1 763 514 4000

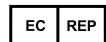

Authorized Representative in the European Community

Medtronic B.V.
Earl Bakkenstraat 10
6422 PJ Heerlen
The Netherlands
+31 45 566 8000

Europe/Middle East/Africa

Medtronic International Trading Sàrl
Route du Molliau 31
Case Postale 84
CH-1131 Tolochenaz
Switzerland
+41 21 802 7000

Australia

Medtronic Australasia Pty Ltd
5 Alma Road
Macquarie Park, NSW 2113
Australia
1800 668 670


Canada

Medtronic of Canada Ltd
99 Hereford Street
Brampton, Ontario L6Y 0R3
Canada
+1 905 460 3800

Technical manuals

www.medtronic.com/manuals

© 2016 Medtronic
M962463A001 B
2016-05-19

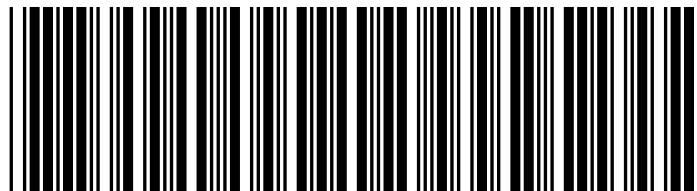

M962463A001