



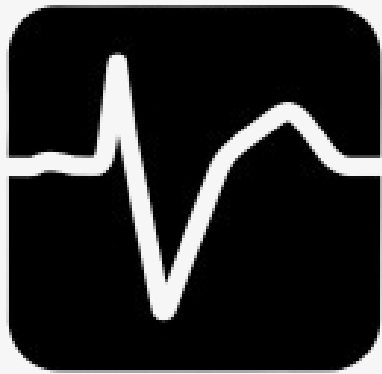

EduPace

User Manual

Quick Start Guide

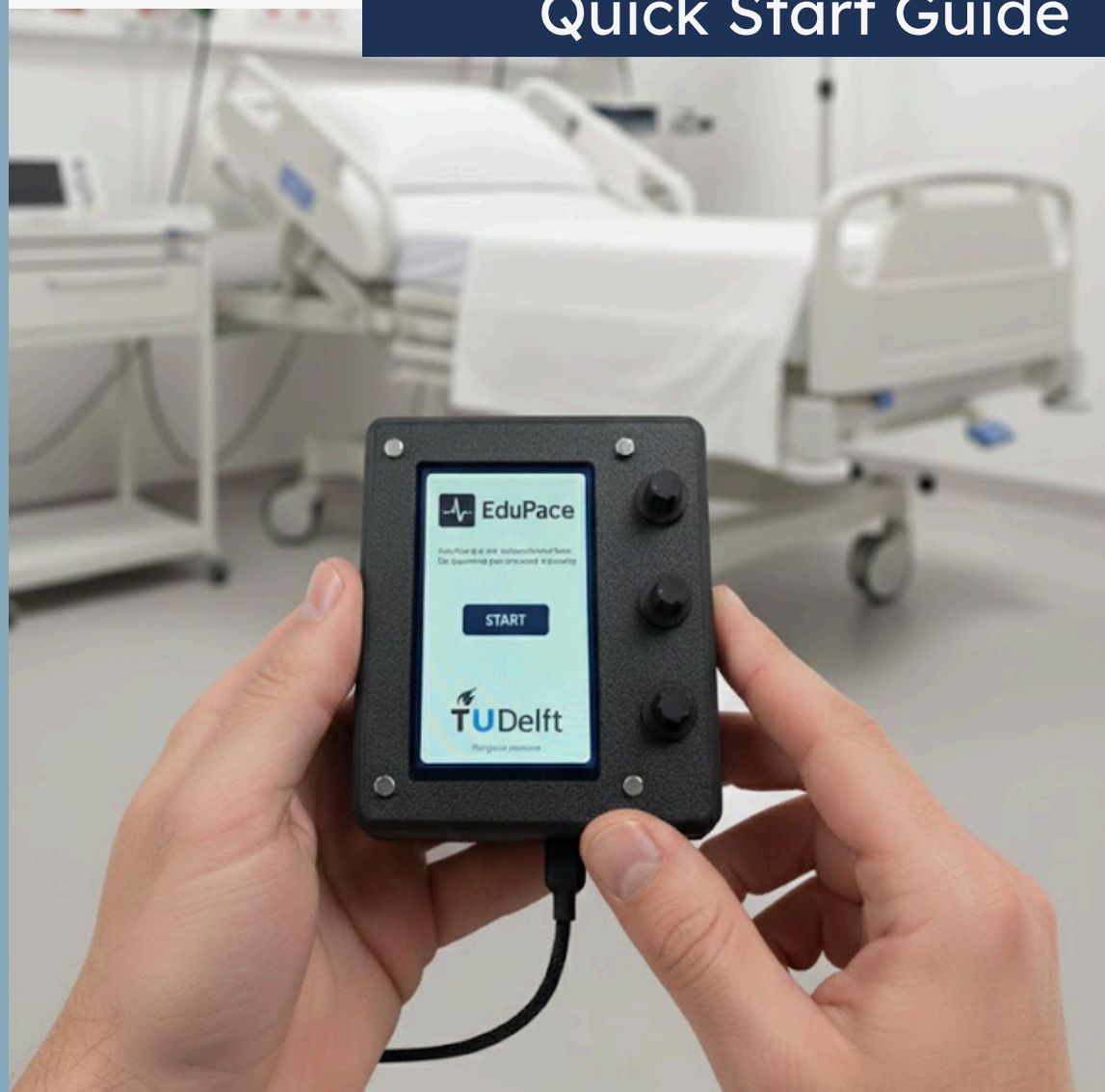



All rights reserved 2025

What is EduPace?

EduPace is an educational training tool designed for practicing the use of an external pacemaker. It has been developed by university students from TU Delft for educational purposes, specifically to support the training of nurses and clinicians.


Reference Device

The design and functionality of EduPace are based on the Medtronic 53401 External Pacemaker, which was used as a reference during development.

Disclaimer

EduPace is not a medical device and is not intended for clinical use. It is an educational exercise tool only and must not be used for diagnosis, treatment, or patient care. Always use certified medical devices in clinical settings.

Controls and Indicators

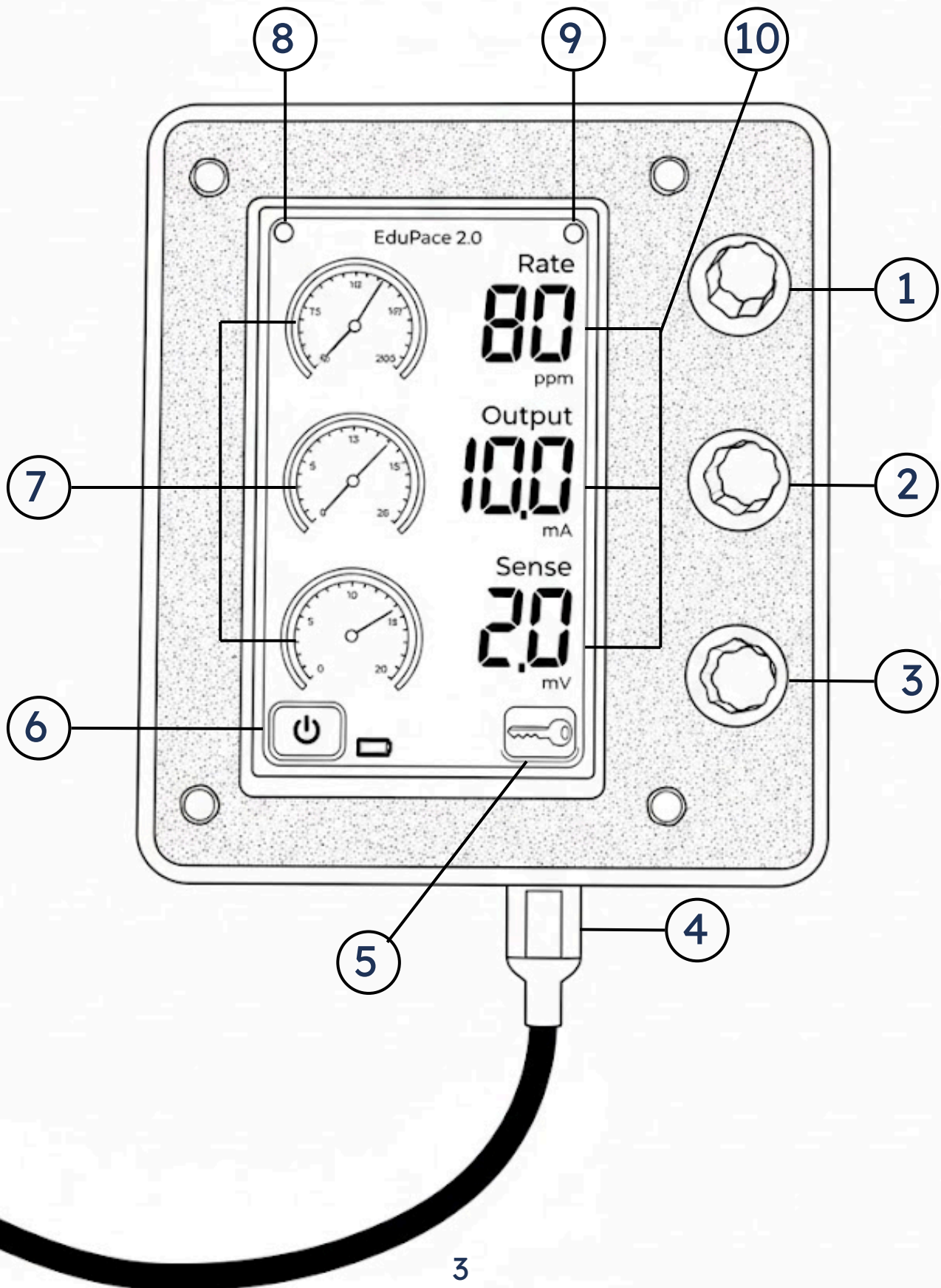


Controls and Indicators

1. RATE dial
2. OUTPUT dial
3. SENSE dial
4. USB-C input
5. Touchable Lock/Unlock button
6. Touchable On/Off button
7. Dashboard for RATE, OUTPUT, SENSE
8. Pacing status indicator LED
9. Sensing status indicator LED
10. Numeric indicators for RATE, OUTPUT, SENSE


System Connector Setup

1. Connect a USB-C cable to the EduPace device and plug the other end into your computer.
2. Launch the EduPace software on your computer and select “Open Training.”

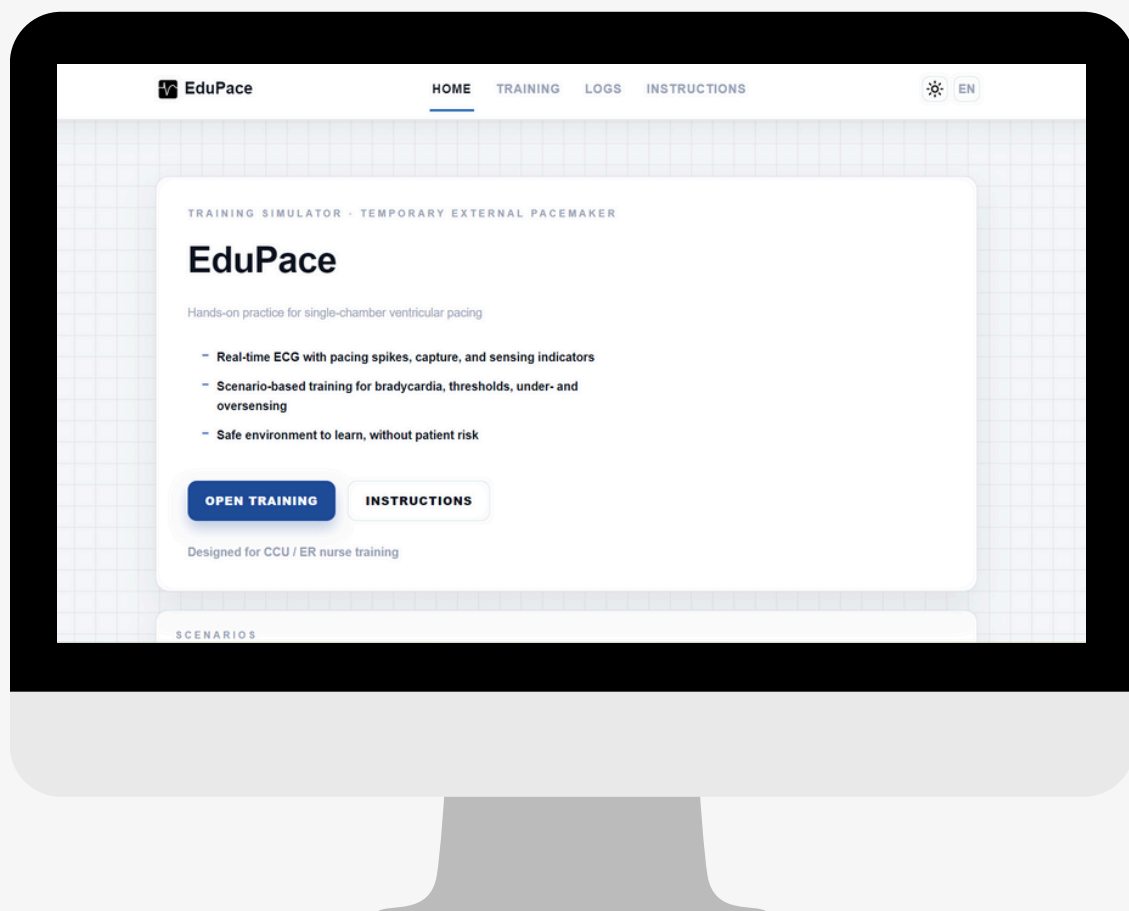

Basic Operation

3. Click “CONNECT” and choose the EduPace device from the list.

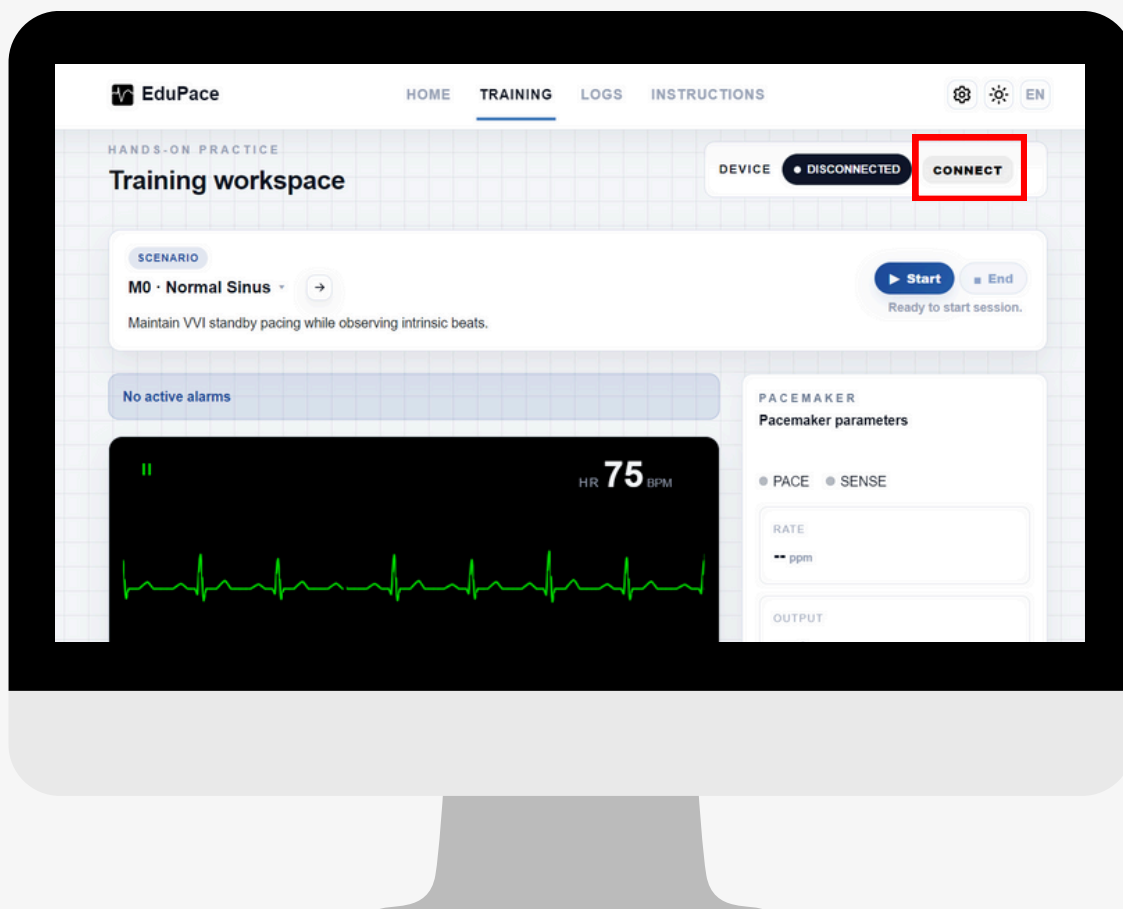

Basic Operation

4. Press START on the EduPace device. The system is now ready for use.



To Turn On/Off

1. To start pacing, briefly press the On/Off button. When the device turns on, the power button lights up green and the background changes from gray to white.



2. To turn off the system, press and hold the On/Off button for more than 2 seconds.

Note: The system must be unlocked before powering off, similar to the Medtronic external pacemaker.

Lock/Unlock

1. To lock the device, briefly press the Lock/Unlock button. The button will turn yellow, indicating that the device is locked. While locked, rotating the RATE, OUTPUT, or SENSE dials will have no effect.



Note: If the RATE, OUTPUT, or SENSE values are not changed for 60 seconds, the system will automatically lock, similar to the Medtronic external pacemaker.

Rate, Output, and Sensitivity Adjustments

1. When EduPace is turned on, the start values are:

RATE: 80 ppm (range 30-200 ppm)

OUTPUT: 10 mA (range 0.1-25 mA)

SENSITIVITY: 2 mV (range 0.1-20 mV)



2. When the system is unlocked, rotate dials to adjust RATE, OUTPUT and SENSE using rotary dials. The left green LED flashes to indicate pacing activity, and the right blue LED flashes to indicate sensing activity.



2025-2026

Contact: EduPace@gmail.com